

PRODUCTS BY TREATMENT SITE



Designed for stereotactic radiosurgery, our products deliver exceptional precision and stability, ensuring accurate targeting for advanced cranial treatments.

Here are our top recommendations for this treatment site. To see more options and variations, check out the full details on the product pages.

BASEPLATES



S-20A
SRS Adaptor
Pg 39

The Klarity SRS Adaptor is designed to fit seamlessly onto a variety of baseplates, making it easy to integrate into your existing setup. Below are some of the compatible baseplates for use with the SRS Adaptor.



R616-SCFA
Uno™ All-In-One System
Pg 30



R630-3HCF
Optek™ S-type Demi-Baseplate for Elekta HexaPOD™
Pg 35



R630-3SCF
Optek™ S-type Overlay
Pg 34

MASKS



S-2000
Klarity Blue™ SRS Open Face Mask Set
Pg 58



PRODUCTS BY TREATMENT SITE



Our breast positioning solutions are tailored for supine treatments, offering exceptional versatility and reliability. Designed for stability and comfort, they ensure precise setups and accommodate a wide range of patient needs.

Here are our top recommendations for this treatment site. To see more options and variations, check out the full details on the product pages.

BASEPLATES



R610-PDCF5
Athena™ Supine Breast System
Pg 40



R616-SCFA
Uno™ All-In-One System
Pg 30



R632-NTCF2
Wingboard
Pg 42

“Thank you for also being so available for us. I know that whenever I have a Klarity question, I will get a response super fast.”

— Kristin Hurocy, RT(T)
Stanford Emanuel Radiation Oncology Center
Turlock, CA

MASKS




R322-5424C1-L
Klarity White™ Head and Thorax Mask
Pg 64



R32-2436
Klarity White™ Athena Breast Mask

PRODUCTS BY TREATMENT SITE

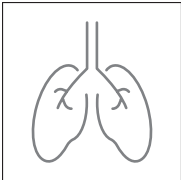

THORAX

Tailored for chest and lung treatments, these systems provide secure and reproducible positioning for improved dosimetric outcomes.

Here are our top recommendations for this treatment site. To see more options and variations, check out the full details on the product pages.

BASEPLATES



R616-SCFA
Uno™ All-In-One System
Pg 30



R610-PDCF5
Athena™ Supine Breast System
Pg 40



R632-NTCF2
Wingboard
Pg 42

Creativity
That Solves
Problems

At Klarity, we approach everything with creativity and a problem-solving mindset. From smarter designs to custom solutions, we think outside the box to help you provide better care with tools that work the way you need them to.

MASKS



RLG461-11ST
Klarity LiteGreen™ S-type
Shoulder Mask
Pg 65



R32-243

PRODUCTS BY TREATMENT SITE

ABDOMEN

Designed to accommodate the unique challenges of abdominal treatments, our solutions deliver optimal patient support and positioning flexibility.

Here are our top recommendations for this treatment site. To see more options and variations, check out the full details on the product pages.

BASEPLATES



R616-SCFA
Uno™ All-In-One System
Pg 30



R632-NTCF2
Wingboard
Pg 42

MASKS



RLG902-2436LA
Klarity LiteGreen™ Uno™
Long Torso Mask
36% perforation
Pg 64

VACUUM BAGS



R7661-20C
58.5cm x 90cm
20-liter fill
Pg 81



R7401-20C
70cm x 70cm
2

PRODUCTS BY TREATMENT SITE



Our pelvic positioning tools provide targeted support and stabilization, offering versatile solutions to enhance accuracy and reproducibility during treatment.

Here are our top recommendations for this treatment site. To see more options and variations, check out the full details on the product pages.

BASEPLATES



R616-SCFA
Uno™ All-In-One System
Pg 30



R634-LCF
Leg Positioner
Pg 43



R634-LCF-2A
Foot Positioner
Pg 44

MASKS



RLG902-2436BA
Klarity LiteGreen™ Uno™ Pelvic Mask
Pg 64



RLG902-3G1
Klarity LiteGreen™ HipFix Mask
Pg 65

KLARITY CUSHIONS™



R550-L2
40cm x 60cm
Pg 72

PRODUCTS BY TREATMENT SITE



BASEPLATES

IN THIS SECTION:

Uno™ All-In-One System.....	30
Optek™ S-type Overlay.....	34
Optek™ S-type Demi-Baseplate for Elekta HexaPOD™	35
Optek™ O-type/U-frame Demi- Baseplate for Elekta HexaPOD™	36
MultiFix™ Baseplate	37
Tilting MultiFix™ Baseplate	38
SRS Adaptor	39
Athena™ Supine Breast System.....	40
Wingboard.....	42
Leg Positioner	43
Foot Positioner	44
Promise™ Proton Therapy System	45



UNO™ ALL-IN-ONE SYSTEM

Precision in Every Position

The Klarity Uno™ is an all-in-one, highly customizable baseplate system suitable for all treatment types.

Available in three sizes, the standard version meets most system requirements, while two specialty versions accommodate narrow bore machines and the Elekta HexaPod system. Each size is also available in an MR-compatible version.

The Uno™ Baseplate set includes a headrest adaptor, right and left handgrips, and a groin lock, and is compatible with standard S-type masks as well as Klarity Uno™ Breast, Thorax, and Pelvic Masks. Separate components are available for different treatment sites with all items offered individually. An SBRT kit, which includes key components for SBRT setup, is also available, along with a storage cart and tool accessory box.

INCLUDED:

- Baseplate
- Headrest Adaptor
- Right and Left Lower Hand Grips
- Groin Lock

Each baseplate is equipped with headrest pins for Silverman headrest compatibility, and sliding bars for Uno™ torso and pelvic mask compatibility. All Uno™ Baseplates have receiving holes for S-type Head and Head and Shoulder Masks.

COMPATIBLE:

S-type Masks	Page 50
Uno™ Breast, Torso, and Pelvic Masks	Page 64
SRS Masks	Page 58
Klarity Cushions™	Page 69
Vacuum Bags	Page 77
Head Supports	Page 91
Indexing Bars	Page 92




UNO™ ALL-IN-ONE SYSTEM



KLARITY UNO


The Uno™ is available in three sizes.

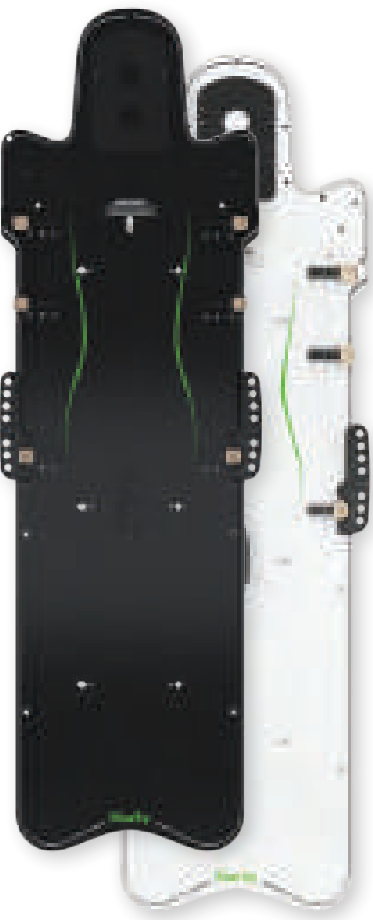

Standard



Narrow

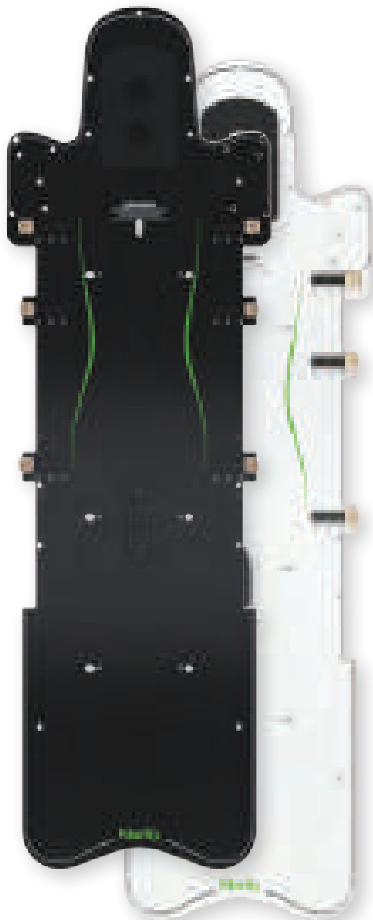

HexaPod

R616-SCFA	R616-SKFA
Carbon Fiber	Kevlar
199.5cm (length) x 61cm (shoulder width) x 54cm (body width)	199.5cm (length) x 61cm (shoulder width) x 54cm (body width)
Max Load 250kg	Max Load 250kg
5.6kg	

R616-NCF	R616-NKF
Carbon Fiber	Kevlar
199.5cm (length) x 54cm (shoulder and body width)	199.5cm (length) x 54cm (shoulder and body width)
Max Load 250kg	Max Load 250kg

R616-SCF1A	R616-SKF1A
Carbon Fiber	Kevlar
199.5cm (length) x 61cm (shoulder width) x 54cm (body width)	199.5cm (length) x 61cm (shoulder width) x 54cm (body width)
Max Load 250kg	Max Load 250kg

COMPATIBLE:

Most Systems

COMPATIBLE:

Small bore CT, MR, and Elekta Unity

COMPATIBLE:

Elekta HexaPod



AVAILABLE COMPONENTS:



R632-NTCF2
Wingboard
Compatible with S-type Masks and the Athena™ Supine Breast System



R616-L1
Knee Positioner
Includes 2 elevation blocks



R616-BE5
Prone Breast Cushion Set
Can be used directly on couch top



R616-WCFA
Overhead Arm Positioner



R616-L3
Foot Positioner



R616-BT1
Uno™ SBRT Set
Pg 33



S20-A
SRS Adaptor
Compatible with Klarity SRS Masks



R616-AIR
Compression Belt



R616-CAB
Uno™ Storage Cart

Brain

Head

SRS

Chest

Lungs

Stomach

SBRT

Hand

Foot

Waves

OPTEK™ S-TYPE OVERLAY

A comfortable, ergonomic S-board for precise head and neck treatments with multiple mask and head support options

The Klarity Optek™ Overlay is a lightweight, durable S-type baseplate designed for head and neck treatments. Constructed from carbon fiber for superior strength and durability, it is also available in a fiberglass version for MR compatibility. Optimized for IMRT and 360° arc treatments, the board offers excellent attenuation and supports both adult and pediatric masks.

Its tapered design enhances patient comfort, while edge indexing ensures precise and reproducible positioning. Adjustable hand grips allow for consistent and comfortable shoulder retraction, further improving treatment accuracy.

- Enhanced indexing
- Extremely lightweight
- Comfortable shoulder retraction
- Strong and durable



R630-3SCF	R630-3SGF
Carbon Fiber	Fiberglass
105cm x 61cm	105cm x 61cm
2.36kg	2.4kg



INCLUDED:

Overlay
Handgrips

ESSENTIAL:

RE-630EA or RE-7E
Klarity Optek™ Indexing Bar

Page 92

COMPATIBLE:

S-type Masks	Page 50
SRS Adaptor	Page 39
SRS Masks	Page 58
Klarity Cushions™	Page 69
Head Supports	Page 91

Brain

Head

SRS

Chest

Lungs

Stomach

SBRT

Hand

Foot

Waves

OPTEK™ S-TYPE DEMI-BASEPLATE FOR ELEKTA HEXAPOD™

Reliable positioning for the head, neck, and shoulders

The Optek™ S-type Demi-Baseplate is tailored for the Elekta HexaPOD™ evo RT System and is also ideal for applications requiring S-type masks with a shortened baseplate. Its ultra-lightweight design features a gently tapered end to enhance patient comfort.

Compatible with all standard S-type masks, it uses standard Silverman pins for precise headrest indexing and supports both adult and pediatric configurations.



R630-3SCF
Carbon Fiber
58cm x 61cm
1.18kg



- Tailored for Elekta HexaPOD™ evo RT System
- Ultra-lightweight design
- Compatible with standard S-type masks

ESSENTIAL:

RE-630EA or RE-7E
Klarity Optek™ Indexing Bar

Page 92

COMPATIBLE:

S-type Masks	Page 50
SRS Adaptor	Page 39
SRS Masks	Page 58
Klarity Cushions™	Page 69
Head Supports	Page 91





OPTEK™ O-TYPE/U-FRAME DEMI-BASEPLATE FOR ELEKTA HEXAPOD™

A lightweight and versatile solution for O-type and U-frame masks

The Klarity Optek™ O-type/U-frame Demi-Baseplate is designed for the Elekta HexaPOD™ evo RT System and is also ideal for applications requiring a shortened baseplate. It is compatible with all O-type 3-point and 5-point masks, as well as U-frame masks.

Constructed from lightweight carbon fiber, the baseplate features a gently tapered lower section to enhance patient comfort.



R630-OHCF
Carbon Fiber
58cm x 61cm
1.26kg

- Compatible with O-type and U-frame masks
- Lightweight carbon fiber construction
- Tapered design for patient comfort

ESSENTIAL: RE-630EA or RE-7E Klarity Optek™ Indexing Bar	Page 92
COMPATIBLE: O-type Masks U-frame Masks Klarity Cushions™ Head Supports	Page 54 Page 56 Page 69 Page 91




MULTIFIX™ BASEPLATE

An easy and flexible solution for the head, neck, and extremities

The Klarity MultiFix™ Baseplate offers a lightweight, durable solution for precise positioning. Constructed from strong carbon fiber, it is fully indexable and compatible with both U-frame and O-type 3-point head masks.

The baseplate supports standard Silverman headrests and includes an optional attachment for prone positioning. In addition to supporting head and neck treatments, the MultiFix™ can be used for extremity positioning.

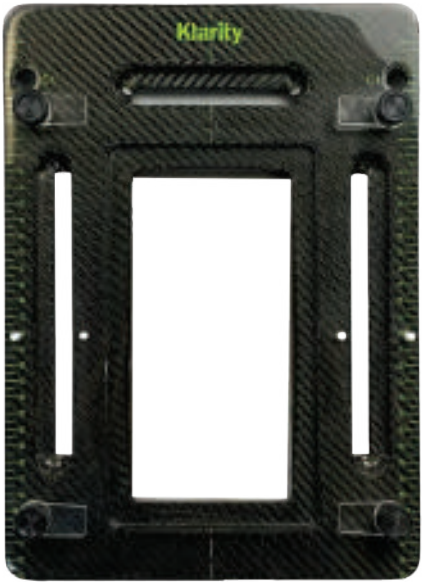

R605-HCF
Carbon Fiber
35.5cm x 25.5cm

- Lightweight and durable
- Compatible with multiple mask types
- Additional indexing along the sides

ACCESSORY:



R605-HPCF
Prone Head MultiFix™
35.8cm x 25.7cm

The Klarity Prone MultiFix™ fits on top of the R605-HCF Baseplate for easy prone positioning. The Prone MultiFix™ is made of lightweight carbon fiber and has cleanable cushions under the forehead and chin.

COMPATIBLE: O-type 3-point Head Masks U-frame Masks Klarity Cushions™ Head Supports	Page 54 Page 56 Page 69 Page 91
--	--



TILTING MULTIFIX™ BASEPLATE

A durable, flexible solution for head tilt positioning

The Klarity Tilting MultiFix™ provides a lightweight, durable solution for precise head tilt positioning, ideal for kyphotic patients. Constructed from strong carbon fiber, it is fully indexable and compatible with O-type and U-frame head masks, as well as standard Silverman headrests.

An optional attachment is available for prone positioning.



- Simple to adjust
- Compatible with multiple mask types
- Additional indexing along the sides



R602-DCF
Carbon Fiber
35.5cm x 25.5cm
250mm maximum height
Tilting Angles: 5°, 10°, 15°, 20°, 25°, 30°, 35°, 40°, 45°

ACCESSORY:



R605-HPCF
Prone Head MultiFix™
35.8cm x 25.7cm

The Klarity Prone MultiFix™ fits on top of the R605-HCF Baseplate for easy prone positioning. The Prone MultiFix™ is made of lightweight carbon fiber and has cleanable cushions under the forehead and chin.

COMPATIBLE:

O-type 3-point Head Masks	Page 54
U-frame Masks	Page 56
Klarity Cushions™	Page 69
Head Supports	Page 91

SRS ADAPTOR

Streamlined Integration for Advanced SRS Positioning



S-20A	S-20B
Plastic	Plastic
30cm x 25cm x 9cm	30cm x 25cm x 9cm
Compatible with Klarity S-type Baseplates	Compatible with other S-type baseplates

- Compatible with all S-type Baseplates
- Easy to attach and detach
- Simple and precise

The Klarity SRS Adaptor is an easy-to-use attachment designed to integrate seamlessly with existing S-type baseplates. It clicks into the standard S-type mask holes of the baseplate, functioning just like a mask.

With the adaptor in place, both top and bottom Klarity SRS masks can be securely attached. Mask height is adjustable from 0–4mm at each pin location using the Simple Shim System™, included with Klarity SRS Masks.

This adaptor is an efficient solution for clinics seeking high-precision SRS treatment without the cost or complexity of a full baseplate system. It is compatible with standard S-type baseplates and Silverman head supports.



COMPATIBLE:

SRS Masks	Page 58
Uno™ All-In-One System	Page 30
Optek™ S-type Overlay	Page 34
Optek™ S-type Demi-Baseplate for Elekta HexaPod™	Page 35
Athena™ Supine Breast System	Page 40
Wingboard	Page 42
All S-type Systems	

ATHENA™ SUPINE BREAST SYSTEM

The ultimate versatility for precise, stable supine breast and torso treatments

The Klarity Athena™ Supine Breast System offers advanced solutions for precise supine breast and torso positioning. The Athena allows for comprehensive arm positioning possibilities with arm cups, upper hand grips, and lower hand grips for asymmetrical positioning. The included Wingboard can be configured with T-bar or U-grip Handles to support customizable shoulder placement. The system also features a securely locking, easily adjustable bum-stop for enhanced stability.

Two Athena™ models are available, each with distinct elevation angles for seamless clinical integration. Both the Athena™ baseplate and Wingboard are compatible with standard S-type head masks, and can also be used with the Klarity SRS Adaptor. Specialized breast and torso masks are also available to meet diverse treatment requirements.



INCLUDED:

- Baseplate
- Wingboard with T-bar and U-grip handles
- Headrest
- Mask Adaptor
- Right and Left Upper Arm Support Cups
- Right and Left Upper Hand Grips
- Right and Left Lower Hand Grips
- Adjustable Bum Stop
- Knee Cushion

COMPATIBLE:

- | | |
|------------------------------|---------|
| S-type Head Masks | Page 50 |
| Athena™ Breast & Torso Masks | Page 64 |
| SRS Adaptor | Page 39 |
| Vacuum Bags | Page 77 |
| Indexing Bars | Page 92 |

R610-PDCF5	R610-PDCF6	R610-PDKF5	R610-PDKF6
Carbon Fiber	Carbon Fiber	Kevlar 	Kevlar 
145cm x 41cm	145cm x 41cm	145cm x 41cm	145cm x 41cm
6.3kg	6.44kg	7.75kg	7.89kg
Max Load 200kg	Max Load 200kg	Max Load 200kg	Max Load 200kg
Tilting Angles: 5°, 10°, 15°, 20°	Tilting Angles: 7.5°, 12.5°, 17.5°, 22.5°	Tilting Angles: 5°, 10°, 15°, 20°	Tilting Angles: 7.5°, 12.5°, 17.5°, 22.5°

ATHENA™ SUPINE BREAST SYSTEM

- Extremely lightweight
- Comprehensive arm placement options
- Compatible with S-type head masks, as well as breast and torso masks



WINGBOARD

A versatile solution with customizable arm positioning and mask compatibility

The Klarity Wingboard is an incredibly versatile baseplate and can be used on its own or indexed directly to both the Uno™ All-In-One System and the Athena™ Supine Breast System. An optional Mask Adaptor allows the use of S-type head masks. It also includes both T-bar and U-grip adjustable handles for customizable arm placement and enhanced patient comfort.



INCLUDED:

- Baseplate
- Mask Adaptor
- T-bar handles
- U-grip handles

COMPATIBLE:

- Uno™ All-In-One System Page 30
- Athena™ Supine Breast System Page 40
- S-type Head Masks Page 50

R632-NTCF2	R632-NTKF2
Carbon Fiber	Kevlar
62cm x 57cm	62cm x 57cm
1.78kg	2.14kg



“We are highly satisfied with the quality and support we get from Klarity. We use Klarity devices across our clinics to ensure the best possible care for our patients.”

— Ainsley Niemkiewicz Fillman, DABR
 Director of Physiscs
 Advocate Radiation Oncology
 Fort Myers, FL

LEG POSITIONER

A comfortable and adjustable solution for precise lower extremity positioning



The Klarity Leg Positioner provides a comfortable, efficient solution for indexing the lower extremities. Designed to attach to any couch top, it integrates seamlessly with other immobilization systems for compatible positioning. Sturdy foam cushions minimize patient fatigue, while the fully adjustable and indexable foot angle ensures precise alignment.

Acrylic hand grips offer comfortable shoulder suppression, with a ball-and-socket locking mechanism allowing reproducible tension in 5cm increments. The system includes two elevation blocks for customizable height adjustment.

Both the Knee Cushion and Elevation Blocks can also be indexed directly to the couch top.

- Indexable baseplate
- Comfortable and intuitive leg positioning
- Reproducible arm straps

R634-LCF	R634-LGF
Carbon Fiber	Kevlar
86.4cm x 44.5cm	86.4cm x 44.5cm
295mm maximum height	295mm maximum height
Tilting Angles: 10°, 20°, 30°	Tilting Angles: 10°, 20°, 30°

INCLUDED:

- Baseplate
- Knee Cushion
- 2 Elevation Blocks
- Foot Holder
- Arm Straps

COMPATIBLE:

- Vacuum Bags Page 77
- Indexing Bars Page 92



FOOT POSITIONER

An adjustable solution for foot positioning with reproducible alignment

The Klarity Foot Positioner offers a comfortable, efficient solution for indexing the feet. It attaches to any couch top and integrates seamlessly with other immobilization systems for compatible positioning.

The system includes a baseplate, foot cushion, and hand grips, with adjustable foot angles ranging from 0–30°. Acrylic hand grips provide comfortable shoulder suppression, while the ball-and-socket locking mechanism ensures reproducible tension in 5cm increments.

- Multiple indexing options
- Easily adjustable angles
- Reproducible arm straps




INCLUDED:

- Baseplate
- Foot Cushion
- Arm Straps

COMPATIBLE:

- Indexing Bars

Page 92

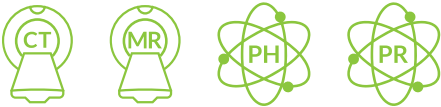

R634-LCF-2A
Acrylic
30cm x 49.4cm
290mm maximum height
Tilting Angles: 10°, 20°, 30°

“As for our work with you and Klarity, you have been nothing but great to us! I would definitely recommend your products to others.”

—McKenzie Gemeiner, B.S., R.T.(T)
US Oncology, Tucson - Craycroft Radiation
Tucson, AZ

PROMISE™ PROTON THERAPY SYSTEM

A precise and lightweight proton therapy solution with user-friendly masks



The Klarity Promise™ Proton System features a lightweight baseplate optimized for proton beam treatments. Adjustable hand grips ensure comfortable, ergonomic shoulder retraction. The system is available as a standard carbon-fiber overlay or an insert-style option for direct couch-top treatments, and is also offered in MR-compatible fiberglass.

Compatible with Klarity's Proton-K Masks, which feature our signature Klarity Green® thermoplastic and a straightforward pin attachment, the Promise™ System provides excellent positioning flexibility. Heat-activated Klarity Cushions™ fit securely in the head cup, offering customized support for the cervical spine.

- Smooth, straightforward mask system
- Comfortable shoulder retraction
- Durable and lightweight

R680-ACF	R680-ACFP1	R680-AGF
Carbon Fiber Overlay	Carbon Fiber Insert	Fiberglass Overlay
120cm x 53cm	157cm x 53cm	120cm x 53cm
2.5kg	4.2kg	2.8kg

INCLUDED:

- Baseplate
- Handgrips

COMPATIBLE:


- Proton-K Masks
- Klarity Cushions™ for Proton Therapy
- Indexing Bars

Page 66
Page 74
Page 92



Streamline Your Workflow

At Klarity, we know time matters. That’s why our products are designed to simplify your workflow and integrate seamlessly with your existing systems. From faster setups to easy adjustments, Klarity tools help you deliver precise, reliable positioning without added complexity. Because when your work flows better, patient care does too.

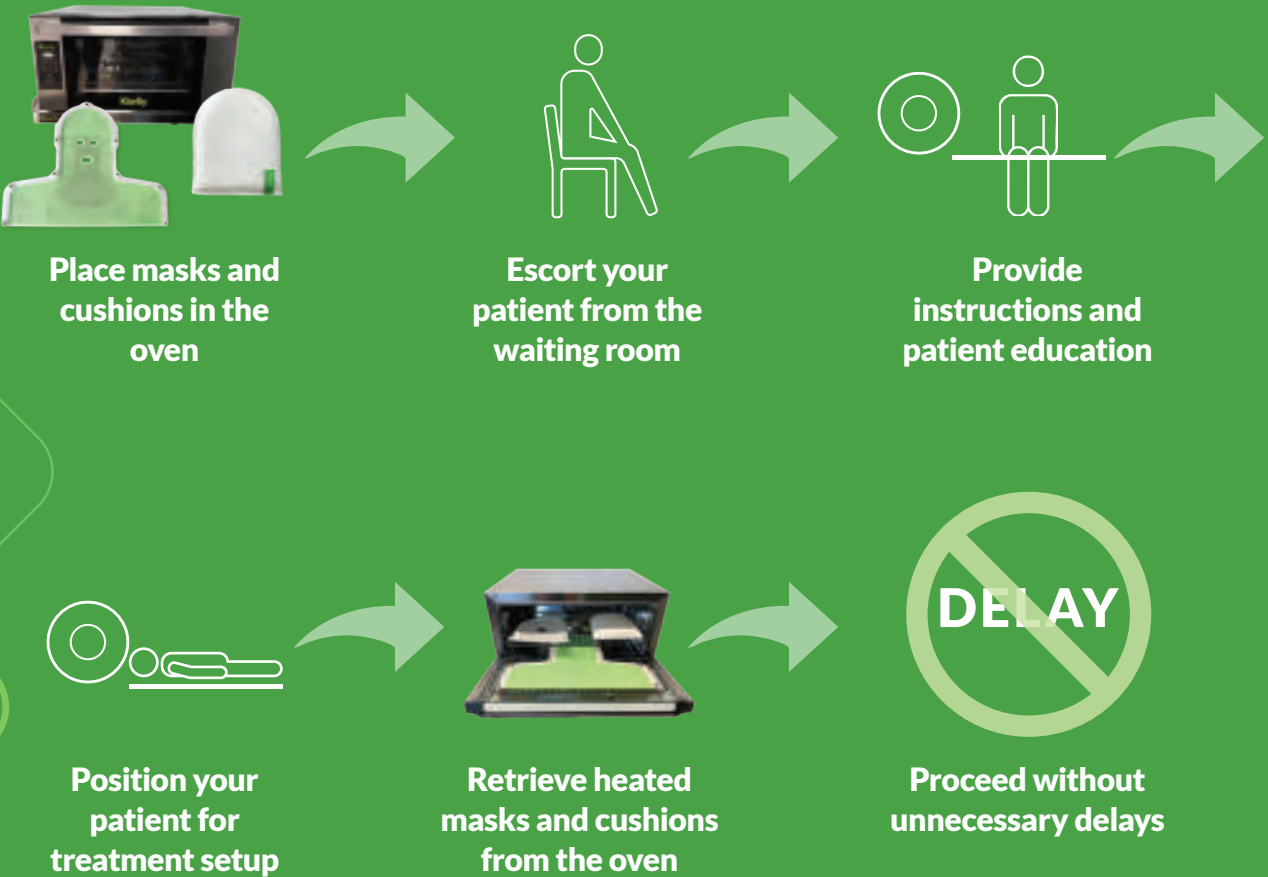

KLARITY THERMOPLASTIC MASKS

IN THIS SECTION:

- Klarity Thermoplastic Masks48
- S-type Head Masks..... 50
- S-type Head and Shoulder Masks 51
- S-type Open Face Head Masks 52
- S-type Open Face Head and Shoulder Masks 53
- O-type Masks 54
- U-frame Masks 56
- SRS Masks..... 58
- V-type Masks 61
- P-type Masks 62
- E-type Masks 63
- Breast, Torso, and Pelvis Masks..64
- Proton Therapy Masks 66

KLARITY THERMOPLASTIC MASKS



Klarity Masks are designed to make every treatment precise, reliable, and efficient. With four advanced thermoplastics—Klarity Blue™, Klarity Green®, Klarity LiteGreen™, and Klarity White™—we offer unmatched strength, reduced shrinkage, and a non-stick coating for easier handling. Each material is optimized for specific frame types and is available in a range of thicknesses, sizes, and perforation patterns to suit diverse clinical needs, ensuring flexibility and accuracy in every setup.

Klarity Masks are easy to prepare, heating at 165°F/73°C in an oven or 160°F/70°C in a water bath. Once warm, the thermoplastic becomes pliable for precise molding and cools to form a rigid, durable fit, delivering consistent performance with every treatment.

Klarity White™

Available for S-type, U-frame, Stereotactic, P-type and Breast, Torso, and Pelvic Masks

Klarity White™ is our original thermoplastic and continues to be the standard for much of the world today. Klarity White™ has a fast set-up for a strong thermoplastic mask. Klarity White™ is less sticky than other standard thermoplastics.

Strong Standard Shrinkage

Molding Time:
45-60 seconds

Cooling Time:
8-10 minutes

Klarity LiteGreen™

Available for O-type, E-type, and Shoulder Masks

Klarity LiteGreen™ is our preferred thermoplastic for 3- and 5-point systems. It has a non-stick coating, so the edges can be pulled apart if they touch, and it drapes nicely for an incredibly easy molding experience. Klarity LiteGreen™ is only available in 2.4mm with a microperforation pattern, which makes this mask stronger than other reinforced thermoplastic.

Stronger Less Shrinkage

Molding Time:
1-2 minutes

Cooling Time:
8-10 minutes

Klarity Green®

Available for S-type, O-type, U-frame, SRS, V-type, P-type, Breast, Torso, and Pelvic Masks, and Proton Therapy Masks

Klarity Green® is our best selling, signature thermoplastic and is our recommended thermoplastic for most applications. It has less shrinkage and has a non-stick coating, so edges can be pulled apart if they connect. Klarity Green® features a slightly longer working-time (30 extra seconds). This means therapists don't have to rush and can create a more carefully molded mask, which allows for a more comfortable patient experience.

Stronger Much Less Shrinkage

Molding Time:
2-3 minutes

Cooling Time:
10-12 minutes

Klarity Blue™

Available for S-type, U-frame, SRS and Proton Therapy Masks

Klarity Blue™ is our strongest thermoplastic, and is ideal for high-precision applications. It also has the least amount of shrinkage. Our blue masks are the most rigid of any mask, with an extremely high resilience once molded and cooled. Klarity Blue™ also sets up quickly, making this ideal for patients who need limited immobilization time, and the blue material drapes gently for smooth, highly accurate mask making.

Strongest Least Shrinkage

Molding Time:
45-60 seconds

Cooling Time:
6-8 minutes



KLARITY MASKS

S-TYPE HEAD MASKS



R461SHT
Klarity White™ S-type Head Mask
42% perforation
3.2mm



RG460SHT
Klarity Green® S-type Head Mask
36% perforation
2.4mm



RB461-3210SHT
Klarity Blue™ S-type Head Mask
10% perforation
3.2mm



R461-4SHT
Klarity White™ S-type Head Mask
AccuPerf™
3.2mm



RG460-4SHT
Klarity Green® S-type Head Mask
AccuPerf™
2.4mm



RG461-5SHT
Klarity Green® S-type Extended Head Mask
AccuPerf™
3.2mm



RB461-3210ST
Klarity Blue™ S-type Head and Shoulder Mask
10% perforation
3.2mm

KLARITY MASKS


S-TYPE HEAD & SHOULDER MASKS



R461ST
Klarity White™ S-type Head and Shoulder Mask
42% perforation
3.2mm



RG460ST
Klarity Green® S-type Head and Shoulder Mask
36% perforation
2.4mm

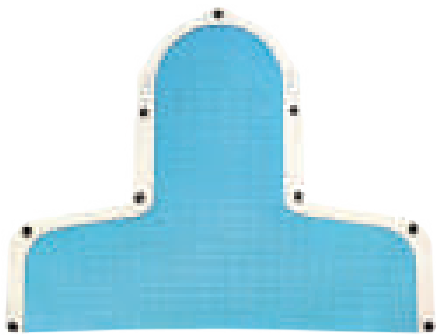

RB461-3210ST
Klarity Blue™ S-type Head and Shoulder Mask
10% perforation
3.2mm



R461-4ST
Klarity White™ S-type Head and Shoulder Mask
AccuPerf™
3.2mm



RG460-4ST
Klarity Green® S-type Head and Shoulder Mask
AccuPerf™
2.4mm



RG460SC
Klarity Green® S-type Pediatric Head and Shoulder Mask
36% perforation
2.4mm
(for R630-3SCF Optek™ S-type Overlay)

RG461-4ST
Klarity Green® S-type Head and Shoulder Mask
AccuPerf™
3.2mm

RG461SC
Klarity Green® S-type Pediatric Head and Shoulder Mask
42% perforation
3.2mm
(for R630-3SCF Optek™ S-type Overlay)

COMPATIBLE:

Uno™ All-In-One System	Page 30
Optek™ S-type Overlay	Page 34
Optek™ S-type Demi-Baseplate for Elekta HexaPod™	Page 35
Athena™ Supine Breast System	Page 40
Wingboard	Page 42
All S-type Systems	



KLARITY MASKS

S-TYPE OPEN FACE HEAD MASKS



R461-8SHTF
 Klarity White™ S-type
 Open Face Head Mask
 AccuPerf™
 3.2mm
 visionrt



RG461-8SHTF
 Klarity Green® S-type
 Open Face Head Mask
 AccuPerf™
 3.2mm
 visionrt



RB461-6SHTF
 Klarity Blue™ S-type
 Open Face Extended Head Mask
 AccuPerf™
 3.2mm
 visionrt



RG461-8SHTM
 Klarity Green® S-type
 Open Eye and Mouth Head Mask
 AccuPerf™
 3.2mm



RB461-6SHTM
 Klarity Blue™ S-type
 Open Eye and Mouth Head Mask
 AccuPerf™
 3.2mm



S-type Open Eye and Mouth Head Mask



S-type Open Face Head Mask

KLARITY MASKS


S-TYPE OPEN FACE HEAD & SHOULDER MASKS



R461-8STF
 Klarity White™ S-type Open Face
 Head and Shoulder Mask
 AccuPerf™
 3.2mm
 visionrt




RG461-8STF
 Klarity Green® S-type
 Open Face Head and Shoulder Mask
 AccuPerf™
 3.2mm
 visionrt

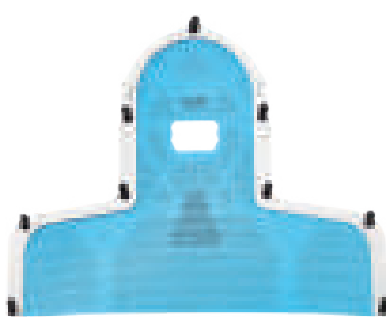

RB461-6STF
 Klarity Blue™ S-type
 Open Face Head and Shoulder Mask
 AccuPerf™
 3.2mm
 visionrt



RG461-8STEM
 Klarity Green® S-type
 Open Eye and Mouth Head and Shoulder Mask
 AccuPerf™
 3.2mm

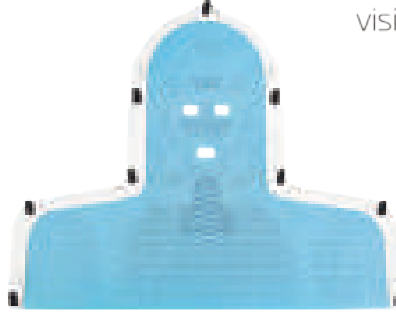

RB461-6STEM
 Klarity Blue™ S-type
 Open Eye and Mouth Head and Shoulder Mask
 AccuPerf™
 3.2mm

Thank you and your team for providing what I believe every center on the Globe should be using. These are the best masks! SPREAD THE WORD!"

—Bob Lemon
 Sutter East Bay Medical Foundation – Sacramento
 Sacramento CA

KLARITY MASKS

O-TYPE 3-POINT HEAD MASKS



RDG105-24A

Klarity Green® O-type 3-point
Head Mask with nosehole
AccuPerf™
2.4mm



RDLG105-2412C

Klarity LiteGreen™ O-type 3-point
Head Mask with nosehole
12% perforation
2.4mm



RDLG105-24AF

Klarity LiteGreen™ O-type 3-point
Open Face Head Mask
AccuPerf™
2.4mm

visionrt



RDG105-32A

Klarity Green® O-type 3-point
Head Mask with nosehole
AccuPerf™
3.2mm



RDLG107-2412C

Klarity LiteGreen™ O-type 3-point
Rounded Head Mask with nosehole
12% perforation
2.4mm



RDLG104-2412C

Klarity LiteGreen™ O-type 3-point
Pediatric Head Mask with nosehole
12% perforation
2.4mm

KLARITY MASKS

O-TYPE 5-POINT HEAD & SHOULDER MASKS



RDG207-24A

Klarity Green® O-type 5-point
Head and Shoulder Mask with nosehole
AccuPerf™
2.4mm



RDLG207-8A24

Klarity LiteGreen™ O-type 5-point
Head and Shoulder Mask with nosehole
12% perforation
2.4mm



RDLG207-24AF

Klarity LiteGreen™ O-type 5-point
Open Face Head & Shoulder Mask
AccuPerf™
2.4mm

visionrt

RDG207-32A

Klarity Green® O-type 5-point
Head and Shoulder Mask with nosehole
AccuPerf™
3.2mm

RDLG207-8A

Klarity LiteGreen™ O-type 5-point
Head and Shoulder Mask with nosehole
12% perforation
3.2mm




RDLG204-2412C

Klarity LiteGreen™ O-type 5-point
Pediatric Head and Shoulder Mask
with nosehole
12% perforation
2.4mm

COMPATIBLE:

- Optek™ O-type/U-frame Demi-Baseplate for Elekta HexaPod™
- Page 36
- MultiFix™ Baseplate
- Page 37
- Tilting MultiFix™ Baseplate
- Page 38
- Orfit Systems

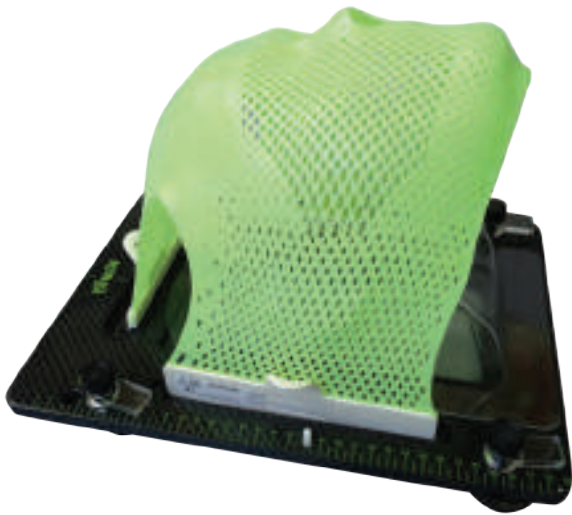

“We love that the Klarity masks are not sticky, and they give us a lot more time to correctly mold them to the patient and to customize them. Overall, it’s a much superior product than what we used in the past.”

—Sarah Steffen, BAS., R.T.(T)
Swedish Cancer Institute – Issaquah
Issaquah, WA



KLARITY MASKS


U-FRAME PROFILEFRAME™ HEAD MASKS



R420W
Klarity White™
ProfileFrame™ Head Mask
36% perforation
2.4mm



RG420W
Klarity Green®
ProfileFrame™ Head Mask
36% perforation
2.4mm

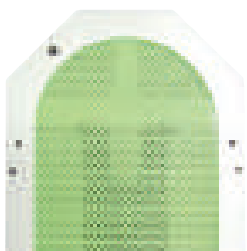

RG420-1W
Klarity Green®
ProfileFrame™ Head Mask
AccuPerf™
2.4mm



RB430W-3210
Klarity Blue™
ProfileFrame™ Head Mask
10% perforation
3.2mm

R430W
Klarity White™
ProfileFrame™ Head Mask
42% perforation
3.2mm

RG430W
Klarity Green®
ProfileFrame™ Head Mask
42% perforation
3.2mm

RG430-1W
Klarity Green®
ProfileFrame™ Head Mask
AccuPerf™
3.2mm



R421W
Klarity White™ ProfileFrame™
Extended Head Mask
36% perforation
2.4mm



RG421W
Klarity Green® ProfileFrame™
Extended Head Mask
36% perforation
2.4mm



RG464-1W
Klarity Green® ProfileFrame™
Extended Head Mask
AccuPerf™
3.2mm



RB431W-3210
Klarity Blue™ ProfileFrame™
Extended Head Mask
10% perforation
3.2mm

R431W
Klarity White™ ProfileFrame™
Extended Head Mask
42% perforation
3.2mm

RG431W
Klarity Green® ProfileFrame™
Extended Head Mask
42% perforation
3.2mm

Also available with
Open Face cutout upon request.
Use ID **RG464-81W**

COMPATIBLE:

Optek™ O-type/U-frame Demi-Baseplate for Elekta HexaPod™	Page 36
MultiFix™ Baseplate	Page 37
Tilting MultiFix™ Baseplate	Page 38
All U-frame Systems	

KLARITY MASKS

U-FRAME HEAD MASKS

EzeFrame™



R430UC
Klarity White™
EzeFrame™ Head Mask
42% perforation
3.2mm

ContourFrame™




R450U
Klarity White™
ContourFrame™ Head Mask
contour perforation
2.4mm

U-Max™ Wide Face Masks



RG420E
Klarity Green®
U-Max™ Wide Face Head Mask
36% perforation
2.4mm

RG430E
Klarity Green®
U-Max™ Wide Face Head Mask
42% perforation
3.2mm

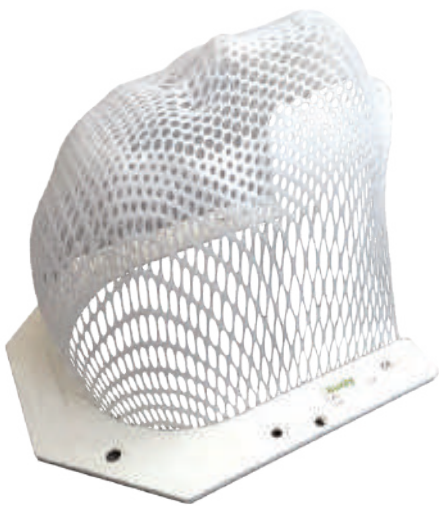

Training Videos at Your Fingertips

Need help? Our training videos are available anytime on our website, with clear, step-by-step instructions to get you up and running fast.

www.klaritymedical.com/instructional-videos




KLARITY MASKS

SRS

Don't get stuck with your SRS Mask!

Klarity SRS Masks provide an extra-strong thermoplastic solution with minimal shrinkage, ensuring precise setups for Stereotactic Radiosurgery.

The two-piece design delivers both anterior and posterior support. Each mask set features removable shims at each pin location, maximizing accuracy and patient comfort with a 0–4 mm range of adjustability. Our premiere Klarity Blue™ thermoplastic offers increased rigidity and stability with the least amount of shrinkage for the most precise treatment alignment.

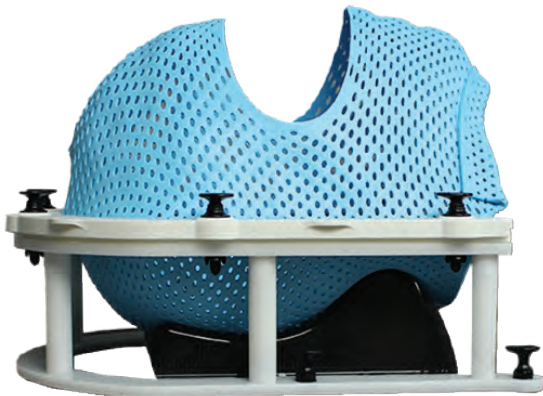

S-2000

Klarity Blue™ SRS Open Face Mask Set

10% perforation
3.2mm

Includes Posterior Mask, Open Face Extended Anterior Mask,
and Simple Shim System™



S-2001

Klarity Blue™ SRS Mask Set

10% perforation
3.2mm

Includes Posterior Mask, Extended Anterior Mask,
and Simple Shim System™

COMPATIBLE:

Klarity SRS Adaptor
Qfix™ Encompass Systems

Page 39




KLARITY MASKS

SRS

Klarity Simple Shim System™

Included with each Klarity SRS Mask

The Simple Shim System™ uses a familiar shimming concept to shim up or down in 1mm increments at each individual pin location. Shims can be added or removed to adjust the height between the posterior and anterior mask from 0mm (no shims) up to 4mm (4 shims). The result is an intuitive shimming technique for secure, customizable treatment.

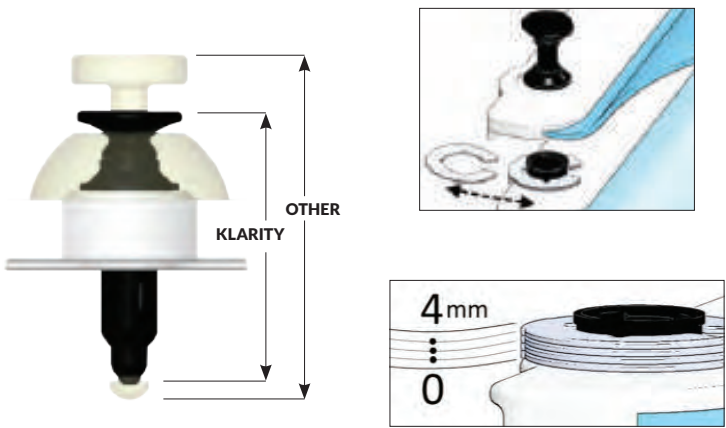

OPTIONAL COMPONENT


The **SRS Single Shim** is placed on the posterior mask and creates an even height across all pin locations in 1mm increments. We recommend starting with 2 shims so you can shim up or down between 0mm-4mm. The Single Shim can be swapped out with the Simple Shim System™.



R-QST-5

SRS Single Shim
1mm



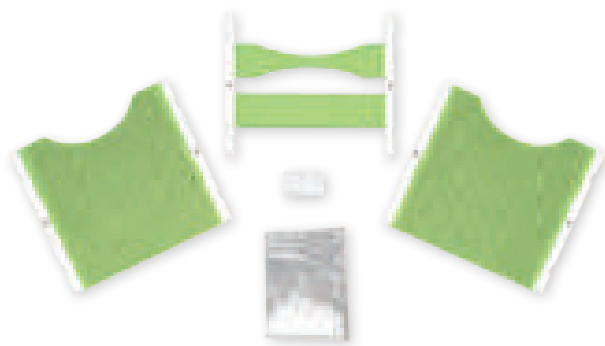

RG408-4AC

Klarly Green® Stereotactic Head Frame Mask Set
Type B

R408

Klarly White™ Stereotactic Head Frame Mask Set
Type B



RG406-1

Klarly Green® Stereotactic Head Frame Mask Set
Type A

R406-1

Klarly White™ Stereotactic Head Frame Mask Set
Type A



RG408-2

Klarly Green® Stereotactic Head and Shoulder Frame Mask Set
Type B

R408-2

Klarly White™ Stereotactic Head and Shoulder Frame Mask Set
Type B
(special order only)

INCLUDED:

Each mask set contains front, middle, and rear masks, a support strip, and a packet of thermoplastic pellets.


COMPATIBLE:

Brainlab Systems



KLARITY MASKS

V-TYPE

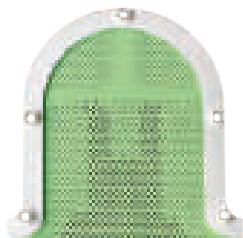

RG460-4VHT

Klarly Green® V-type Head Mask
AccuPerf™
2.4mm



RG460-4VT

Klarly Green® V-type Head and Shoulder Mask
AccuPerf™
2.4mm

COMPATIBLE:

Bionix Versa Systems



“Loving the masks! It is so nice for my patients not to get their hair pulled by the mask or worse, beards. Klarly has amazing products to offer.”

— Renae Kozak, BSRT(R)(T)
Northwestern Medicine Lake Forest Hospital Cancer Center
Lake Forest, IL

KLARITY MASKS

P-TYPE



RG-PRT3-24A

Klarity Green® P-type 3-point Head Mask with nosehole

AccuPerf™

2.4mm



RG-PRT5-24A

Klarity Green® P-type 5-point Head and Shoulder Mask with nosehole

AccuPerf™

2.4mm



R-PRT2-2415

Klarity White™ P-type 2-point Head Mask with nosehole

15% perforation

2.4mm



R-PRT33-2415

Klarity White™ P-type 3-point Short Head Mask with nosehole

15% perforation

2.4mm



R-PRT3-2415

Klarity White™ P-type 3-point Head Mask with nosehole

15% perforation

2.4mm



R-PRT31-2415

Klarity White™ P-type 3-point Rounded Head Mask with nosehole

15% perforation

2.4mm



R-PRT44-2415

Klarity White™ P-type 3-point Extended Head Mask with nosehole

15% perforation

2.4mm



R-PRTL-2415

Klarity White™ P-type 4-point Head and Neck Mask with nosehole

15% perforation

2.4mm



R-PRT55-2415

Klarity White™ P-type 5-point Narrow Head and Shoulder Mask with nosehole

15% perforation

2.4mm



R-PRT5-2415

Klarity White™ P-type 5-point Head and Shoulder Mask with nosehole

15% perforation

2.4mm

COMPATIBLE:

CQ Posifix® Systems



KLARITY MASKS

E-TYPE MASKS



RLG-ESH3-2412C

Klarity LiteGreen™ E-type Esprit and Icon™ Head Mask with nosehole

12% perforation

2.4mm



RCLG-ESH3-2412CMB

Klarity LiteGreen™ E-type Esprit and Icon™ Head Mask with 25mm nosehole

12% perforation

2.4mm




RCLG-ESH3-2412CN

Klarity LiteGreen™ E-type Esprit and Icon™ Head Mask with 13mm triangle nosehole

12% perforation

2.4mm

Masks for Elekta Fraxion™

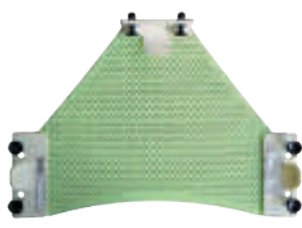

RLG-ESH1-2412C *

Klarity LiteGreen™ E-type Fraxion™ Head Mask

12% perforation

2.4mm



RLG-ESH2-2412C *

Klarity LiteGreen™ E-type Fraxion™ Rounded Head Mask

12% perforation

2.4mm

Masks for Elekta HeadSTEP™



RLG-ELH-2412CN *

Klarity LiteGreen™ E-type HeadSTEP™ Head Mask

22% perforation

2.0mm



RLG-ELS-2412CN *

Klarity LiteGreen™ E-type HeadSTEP™ Head and Shoulder Mask

22% perforation

2.0mm

COMPATIBLE:

Klarity Cushions™ for Elekta Esprit and Icon™ Systems

Elekta's Esprit, Icon™, Fraxion™, and HeadSTEP™ Systems





KLARITY MASKS

BREAST, TORSO & PELVIC MASKS

Uno™ Masks



RLG322-S42436

Klarity LiteGreen™ Uno™ Head and Thorax Mask
36% perforation
2.4mm
59cm x 46cm

R322-S42436

Klarity White™ Uno™ Head and Thorax Mask
36% perforation
2.4mm
59cm x 46cm



RLG902-2436LA

Klarity LiteGreen™ Uno™ Long Torso Mask
36% perforation
2.4mm
59cm x 46cm

R902-2436LA

Klarity White™ Uno™ Long Torso Mask
36% perforation
2.4mm
59cm x 46cm



RLG322-S424C-R

Klarity LiteGreen™ Uno™ Head and Breast Mask
Right Side
36% perforation
2.4mm
59cm x 46cm

R322-S424C-R

Klarity White™ Uno™ Head and Breast Mask
Right Side
36% perforation
2.4mm
59cm x 46cm



RLG902-2436TL

Klarity LiteGreen™ Uno™ Thorax Mask
36% perforation
2.4mm
41cm x 46cm

R902-2436TL

Klarity White™ Uno™ Thorax Mask
36% perforation
2.4mm
41cm x 46cm



RLG322-S424C-L

Klarity LiteGreen™ Uno™ Head and Breast Mask
Left Side
36% perforation
2.4mm
59cm x 46cm

R322-S424C-L

Klarity White™ Uno™ Head and Breast Mask
Left Side
36% perforation
2.4mm
59cm x 46cm



RLG902-2436BA

Klarity LiteGreen™ Uno™ Pelvic Mask
36% perforation
2.4mm
56cm x 46cm

R902-2436BA

Klarity White™ Uno™ Pelvic Mask
36% perforation
2.4mm
56cm x 46cm

COMPATIBLE:

Uno™ All-In-One System Page 30

Athena™ Masks



R322-S424C1-L

Klarity White™ Athena™ Head & Thorax Mask
Left Side
3.2mm



R322-S424C1-R

Klarity White™ Athena™ Head & Thorax Mask
Right Side
3.2mm



R32-2436

Klarity White™ Athena™ Breast Mask
2.4mm

COMPATIBLE:

Athena™ Supine Breast System Page 40

KLARITY MASKS

BREAST, TORSO & PELVIC MASKS

Shoulder Mask



RLG461-11ST

Klarity LiteGreen™ S-type Shoulder Mask AccuPerf™
3.2mm

COMPATIBLE:

Uno™ All-In-One System Page 30
Optek™ S-type Overlay Page 34
Optek™ S-type Demi-Baseplate for Elekta HexaPod™ Page 35
All S-type Systems

O-type Breast & Pelvic Masks

Special Order Only



RD305-3242C

Klarity White™ O-type 2-point Breast Mask
42% perforation
3.2mm
25cm x 40.8cm

RD301-3242

Klarity White™ O-type 2-point Breast Mask
42% perforation
3.2mm
34cm x 23cm



RD901-2436

Klarity White™ O-type 4-point Pelvic Mask
36% perforation
2.4mm
46.3cm x 48cm

RD901-3242

Klarity White™ O-type 4-point Pelvic Mask
42% perforation
3.2mm
46.3cm x 48cm

COMPATIBLE:

Orfit Systems

PelvicFix™ Masks



R902-3A

Klarity White™ PelvicFix™ Mask
42% perforation
3.2mm
46cm x 56cm



RLG902-3G1

Klarity LiteGreen™ HipFix Mask
Reinforced perforation
3.2mm
58.8 x 45.7cm

COMPATIBLE:

Klarity Classic Series: R606-ACD Klarity PelvicFix™ System
Civco HipFix System

Breast Masks

Special Order Only



R91-2436M

Klarity White™ Pelvic Mask
36% perforation,
2.4mm
46cm x 48cm

R91-3242M

Klarity White™ Pelvic Mask
42% perforation
3.2mm
46cm x 48cm

R91-2400M

Klarity White™ Pelvic Mask Solid
2.4mm
46cm x 48cm

COMPATIBLE:

Klarity Classic Series: R611-5SDCF Supine Breast System



KLARITY MASKS PROTON THERAPY

Proton-K Masks



RG-PH-2436C

Klarly Green® Proton-K
Head Mask
36% perforation
2.4mm



RG-PS-2436C

Klarly Green® Proton-K
Head and Shoulder Mask
12% perforation
2.4mm

RG-PH-3242C

Klarly Green® Proton-K
Head Mask
42% perforation
3.2mm

RG-PS-3242C

Klarly Green® Proton-K
Head and Shoulder Mask
42% perforation
3.2mm

COMPATIBLE

Klarly Promise™ Proton System Page 45
Klarly Cushions™ for Proton Therapy Page 74

“The variety of masks
offered covers a multitude
of treatment modalities
without sacrificing quality,
reproducibility or comfort.”

—Ruth Barbir, R.T.(T),
South Florida Proton Therapy Institute
Delray, FL

Proton-B Masks



RB-BSL-3210

Klarly Blue™ Proton-B
Large Head and Shoulder Mask
10% perforation
3.2mm
48cm width



RG-BSL-3242

Klarly Green® Proton-B
Large Head and Shoulder Mask
42% perforation
3.2mm
48cm width



RG-BHM-3242

Klarly Green® Proton-B
Medium Head Mask
42% perforation
3.2mm
31cm width



RG-BHS-3242

Klarly Green® Proton-B
Small Head Mask
42% perforation
3.2mm
26cm width



RG-BHS-A2436

Klarly Green® Proton-B
Pediatric Head Mask
36% perforation
2.4mm
24cm width

COMPATIBLE

Qfix BoS™ System



5 Reasons Clinics Love Klarity Cushions™

- 1. No Odor** – Clean and odor-free for a better treatment experience.
- 2. No Leaks** – No liquid components, no risk of leaks—ever.
- 3. No Waste** – Heat-activated only, so they’re safe from punctures or torn packaging.
- 4. Easy to Modify** – Simple to adjust, reshape, or refine as needed.
- 5. Customizable** – Easily molded to provide personalized patient support.

KLARITY CUSHIONS™

IN THIS SECTION:

Klarity Cushions.....	70
Sizes and Specifications	71
Specialty Cushions:	
Klarity Cushions™ for Proton Therapy.....	74
Klarity Cushions™ for Elekta Esprit & Icon™ Systems	74
Tapered Cushions	75
Cushion Accessories	75

KLARITY CUSHIONS™

Odorless, heat-activated customized support

Klarity Cushions™ provide secure, personalized support for the head, neck, and extremities, making them ideal for cranio-spinal, head and neck, and custom positioning needs. These low-density thermoplastic molds are designed to soften when heated and become firm upon cooling, ensuring precise placement for accurate treatment.

Odorless and dry, Klarity Cushions™ are easy to prepare and offer a working time of 2–4 minutes, fully setting in 5–7 minutes. They are patient-specific but reheatable and reformable, providing exceptional flexibility for therapists. With low attenuation and unmatched comfort, Klarity Cushions™ are a trusted choice worldwide.

Available in various shapes, sizes, and thicknesses, Klarity Cushions™ support a wide range of positioning needs, including head, neck, shoulders, torso, and extremities. Custom designs are also available upon request to meet unique clinical requirements.

- Odorless
- Heat-Activated
- Extremely firm and reproducible
- Remoldable



KLARITY CUSHIONS™ U-CUSHIONS

 R550-U 23cm x 30cm	 R550-U2 23cm x 30cm with additional fill	 R550-U3 23cm x 35cm
 R550-U4 20cm x 25cm	 R550-U5 23cm x 25cm	 R550-U6 23cm x 38cm with additional fill
 R550-U7 23cm x 35cm with double fill	 R550-U8 23cm x 38cm with triple fill	 R550-U9 23cm x 38cm



KLARITY CUSHIONS™

KLARITY CUSHIONS™



R550-S
15cm x 20cm



R550-M
20cm x 25cm



R550-M1
23cm x 30cm



R550-M2
20cm x 25cm
with double fill



R550-M3
23cm x 35cm



R550-L
20cm x 40cm



R550-L2
40cm x 60cm

“Klarity cushions are the best in the game! They can be reshaped if necessary and the extra fill provides the comfort our patients want and need!”

—**Ruth Barbir, R.T.(T)**
South Florida Proton Therapy Institute
Delray, FL



R550-T
46cm x 41cm



R550-TL
46cm x 46cm



R550-T2
46cm x 41cm
with additional fill



R550-T2L
46cm x 46cm
with additional fill



R550-H1
26cm x 35cm



R550-HS2L
46cm x 50cm
with additional fill

Pediatric
Klarity Cushion™



R550-P
28cm x 36cm

KLARITY CUSHIONS™ SPECIALTY

Klarity Cushions™ for Proton Therapy



R550-BOS
46cm x 45cm
with additional fill



R550-BOS3
48cm x 63cm
with additional fill



R550-K
44cm x 44cm
with tapered fill



R550-CM1
15cm x 25cm
with additional fill



R550-PS
16cm x 31cm

Klarity Cushions™ for Elekta Esprit & Icon™ Systems



R550-CEL
for Elekta Esprit and Icon™ Systems
17.5cm x 35cm




R550-CELC
for Elekta Esprit and Icon™ Systems
17.5cm x 35cm
1/2 fill level



R550-CELC1
for Elekta Esprit and Icon™ Systems
17.5cm x 35cm
1/3 fill level

KLARITY CUSHIONS™ ACCESSORIES AND TOOLS

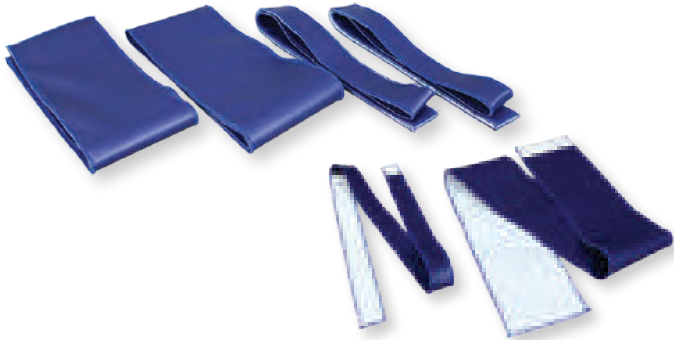

R565-605C

Klarity Cushion™ Cleanable Strap Set
for Cushion Molding

The cleanable strap set is a positioning aid that makes it easier for a single therapist to mold a Klarity Cushion™. The straps are made from a comfortable and stretchable material that is easy to clean. Velcro strips attach to the Klarity Cushion™ and sides of the couch, providing comfortable but firm suppression during setup. Straps are used for setup only and should be removed before treatment.

- Each Strap Set comes with the following:
 - Two 5cm x 152cm straps
 - Two 15cm x 152cm straps
 - Two adhesive tabs for the couch



R560-S

Klarity Shell Mold for S-type and U-frame

R560-O

Klarity Shell Mold for O-type

R560-M

Klarity Shell Mold for R605-HCF MultiFix™ Baseplate
The Klarity Shell Mold is a reusable frame for molding the H1 and HS2L Klarity Cushions™. Once the cushion cools within the frame, the Shell Mold should be removed before using a mask. Includes three head sizes for each frame type.

Tapered Cushions



R550-W1

31cm x 23cm
9cm max height
with tapered fill

Custom Products Made Easy

From custom mask cutouts to unique vacuum bag and cushion sizes, Klarity makes custom work simple, affordable, and accessible. No stress. Just the right solution for your clinic’s unique needs.

If you need something specific, we’re ready to make it happen. Just ask.

sales@klaritymedical.com

KLARITY VACUUM BAGS

IN THIS SECTION:

Vacuum Bags 78

Vacuum Bag Accessories 79

Sizes and Specifications 80

KLARITY VACUUM BAGS

A stronger, more durable vacuum bag



Klarity Vacuum Bags are designed for strength, reliability, and long-term use. A durable, non-permeable nylon shell prevents punctures and tears, ensuring each bag is completely air-tight and leak-proof. Filled with low-density polystyrene mini-spheres, they provide secure, consistent support.

Easy to clean, Klarity Vacuum Bags are built for multi-patient use and require less frequent replacement thanks to their durable construction. Replaceable valves and adaptable connectors ensure compatibility with any vacuum pump, while a corner grommet allows for convenient storage.

Adaptors

Klarity Vacuum Bags are automatically equipped with valves for use with Klarity or Civco Vacuum Pumps. Please inquire if additional Vacuum Pump adaptors are needed.

R7800-QAD Adaptor for QFix™ Vacuum Pump <i>Allows QFix™ Pump to be used with Klarity Vacuum Bags</i>	R7800-QADF Adaptor for Klarity Vacuum Pump <i>Allows Klarity Pump to be used with QFix™ Vacuum Bags</i>	R7800-EAD Adaptor for Elekta Vacuum Pump <i>Allows Elekta Pump to be used with Klarity Vacuum Bags</i>	R7800-NAD Adaptor for Bionix™ Vacuum Pump <i>Allows Bionix™ Pump to be used with Klarity Vacuum Bags</i>	R7800-NADF Adaptor for Klarity Vacuum Pump <i>Allows Klarity Pump to be used with Bionix™ Vacuum Bags</i>
---	--	---	---	--

KLARITY VACUUM BAGS ACCESSORIES



KSC-3
Klarity Storage Cart
45.7cm x 121.9cm x 182.8cm
The cart holds up to 227 kg (500 lbs), features locking wheels, and includes three adjustable shelves. Movable rungs accommodate vacuum bag grommets for easy storage.



KVP-500
Klarity Vacuum Pump with AirSwitch™
220v or 110v
Includes AirSwitch™, 2 valves, and 2 self-retracting vacuum tubes that are 3m (10ft) in length.



KVP-100
Klarity AirSwitch™
The AirSwitch™ is a lightweight handheld control with three settings—In, Hold, and Out—for precise bag adjustments. It connects to your existing vacuum pump and works with all vacuum bag types and valves. Please specify your pump type when ordering.



R7800-VBI
Vacuum Bag Indexing Bar Set
Our Vacuum Bag Indexing Bars allows you to transform any vacuum bag into an indexed bag. The Vacuum Bag Indexing Bars attach to your existing indexing bar and features 4 positions so you can easily adjust the indexing width.



KLARITY VACUUM BAGS

Klarity Vacuum Bags are available in the following sizes and specifications.
All bags can be ordered with an indexing bar. Please specify valve type when requesting a quote. Inquire directly for custom sizes.



R7100-01C
25cm x 20cm
1-liter fill



R7101-03C
40cm x 30cm
3-liter fill



R7106-2.25C
50cm x 25cm
2.25-liter fill



R7301-15C
70cm x 50cm
15-liter fill



R7401-20C
70cm x 70cm
20-liter fill



R7506-40C
100cm x 100cm
40-liter fill



R7504-17.5C
100cm x 70cm
17.5-liter fill



R7504-30C
100cm x 70cm
30-liter fill



R7504-35C
100cm x 70cm
35-liter fill



R7504-35BC
100cm x 70cm
35-liter fill
with Indexing Bar



R7505-30C
100cm x 80cm
30-liter fill



R7536-30C
120cm x 70cm
30-liter fill

KLARITY VACUUM BAGS



R7540-30C
150cm x 75cm
30-liter fill



R7540-60C
150cm x 75cm
60-liter fill



R7512-55C
150cm x 80cm
55-liter fill



R7509-75C
150cm x 100cm
75-liter fill



R7200-05C
50cm x 70cm
5-liter fill



R7202-05C
67cm x 62cm
5-liter fill



R7626-20C
53cm x 104cm
20-liter fill



R7661-20C
58.5cm x 90cm
20-liter fill



R7644-32.5BC
90cm x 36.5cm
32.5-liter fill
with Indexing Bar



R7131-5NLTC
SBRT Knee Bridge Vacuum Bag
28cm x 80cm
5-liter fill



R7524-120C
200cm x 100cm
120-liter fill

“During my 20 years of experience in Radiation Oncology, Klarity has become the leader in immobilization devices. Klarity provides products and clinical support that are key to improving patient outcomes and increasing the opportunities to cure cancer. Not only does Klarity make exceptional products, they provide educational opportunities for Radiation Oncologists, Medical Dosimetrists and Radiation Therapists all over the world.”

—Chuck Lemmond, B.S. RT(R)(T)(ARRT)
College of Veterinary Medicine- Auburn University
Auburn, AL

ACCESSORIES

IN THIS SECTION:

The Klarity BiteLok®	84
Bolus	86
Patient Specialized Accessories:	
Shields	88
Skin Marking.....	89
General Treatment Accessories:	
Compression Belt.....	90
Head Supports	91
Indexing Bars	92
Table Pads & Wedges.....	94
Straps.....	95

THE KLARITY BITELOK®

PRECISION MEETS COMFORT




- Breathable, low-profile mouthpiece
- Secure indexing to mask
- Removable for easy cleaning
- Optional tongue depressor/diverter
- Comfortable for smaller mouths

The Klarity BiteLok® is a lightweight, breathable mouthpiece designed for secure and reproducible head and neck positioning. Using fast-setting dental putty—no water or warming required—the BiteLok quickly indexes to the patient’s mouth. A hands-free Quickfit™ tool molds a thermoplastic mask around the edges of the BiteLok, creating a precise connection between the patient and the mask. The mouthpiece can be easily removed for cleaning.

With its low-profile design and wide airway passage, the BiteLok offers superior comfort, even for patients with smaller mouths. It’s an innovative solution for accurate, reliable positioning in head and neck treatments, delivering both precision and patient comfort.



QuickFit tool indexing a thermoplastic mask to the BiteLok

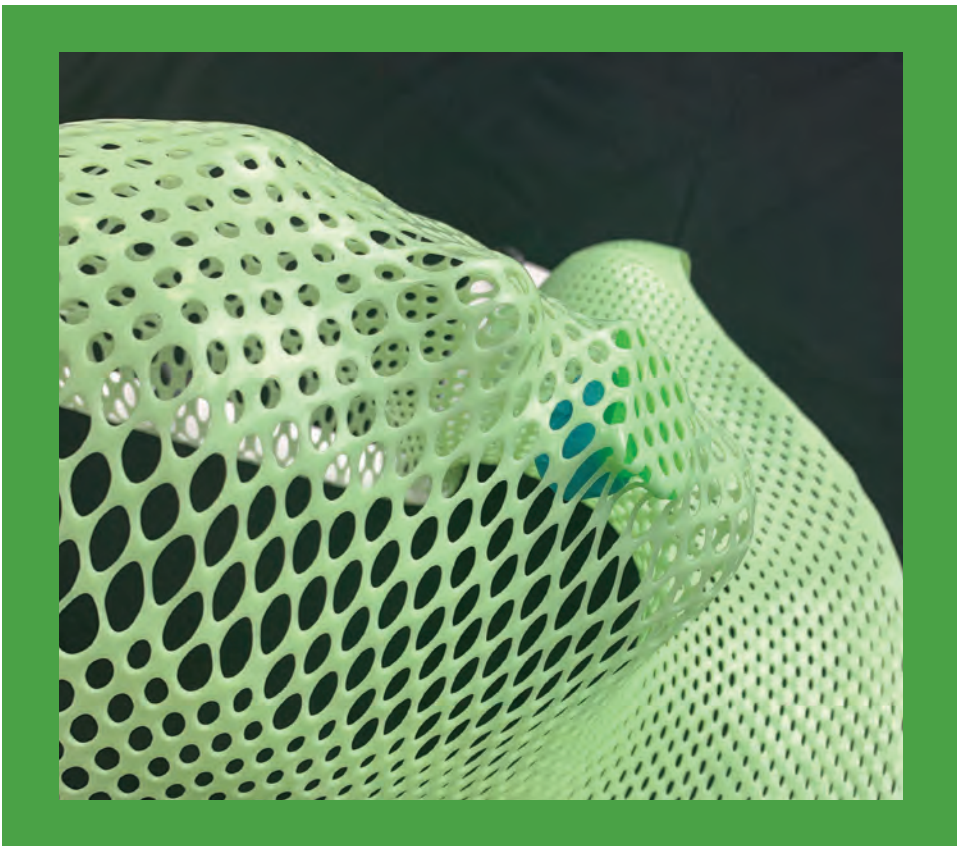

*US and International patents pending



THE KLARITY BITELOK®

BiteLok® Kit



RNB-100
BiteLok Kit

The BiteLok Kit comes as a complete package with all tongue placement options included and contains:

- BiteLok
- Two Putty Packets
- QuickFit Tool
- Short and Long Tongue Depressors
- Right and Left Tongue Diverters
- Tongue Device Removal Tool

BiteLok®



RNB-10
BiteLok

The BiteLok can be ordered as a simplified package containing:

- BiteLok
- One Putty Packets
- QuickFit Tool



Tongue Placement Kit

Optional tongue depressors and diverters provide added stability and customization.

RNB-11 Tongue Placement Accessory Kit

The Tongue Placement Accessory Kit also be ordered separately and contains:

- Short and Long Tongue Depressors
- Right and Left Tongue Diverters
- Tongue Device Removal Tool

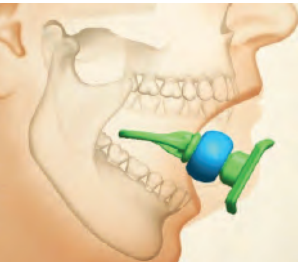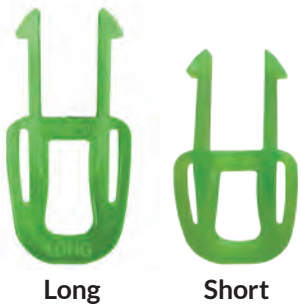

The BiteLok offers optional tongue depressors and diverters for added stabilization and mouth placement options. These devices snap securely into place and are easily removed using the included removal tool.



BOLUS

Klarity offers four types of standard bolus:

Klarity Bolus™: a filmless, gel-based bolus with the optimal amount of stickiness	Bolx™ Bolus: a tissue-equivalent bolus encased in a layer of polyethylene skin	Bolus Pellets: heat-activated pellets ideal for irregular contours, such as the nose or ear	Klarity ThermoSheets™: easy to form, heat-activated thermoplastic sheets for bolus applications
---	--	---	---

Klarity Bolus™

Klarity Bolus™ is a filmless, tissue-equivalent bolus designed for radiation therapy. Its flexibility allows it to conform closely to the patient's anatomy, and it has just the right amount of tackiness for a secure fit. The bolus is translucent, enabling clear visualization of portal marks, treatment areas or lesions underneath.


Klarity Bolus™ can be cut or layered to achieve the desired thickness, and is constructed from silicone gel reinforced with a thin polyester mesh, providing both dimensional stability and reliable performance in every application.



RKB-3003 30cm x 30cm 3mm	RKB-4005 40cm x 40cm 5mm
RKB-3005 30cm x 30cm 5mm	RKB-4010 40cm x 40cm 10mm
RKB-3010 30cm x 30cm 10mm	

Bolx™ Bolus

Bolx I™ is tissue-equivalent flexible build up materials. It is designed to increase the surface dose and offers excellent conformity to body contours. Bolx I™ is made of a cured transparent gel and is encased in a layer of polyurethane skin.

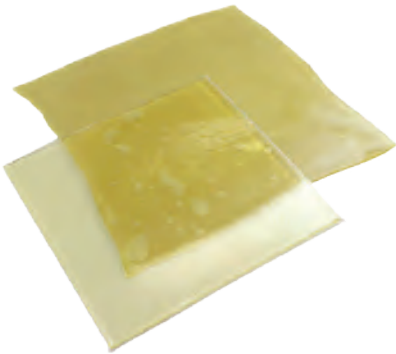

KAC-31033030 30in x 30in 3mm	KAC-31102 30in x 30in 10mm
KAC-31034040 40in x 40in 3mm	KAC-31152 30in x 30in 15mm
KAC-31052 30in x 30in 5mm	KAC-31202 30in x 30in 20mm

BOLUS

Bolus Pellets

Klarity Bolus Pellets are an ideal choice for custom bolus applications. Once heated in an oven or water bath, they become transparent and easy to mold into any desired shape—particularly useful for areas like the nose and ear. Their dosimetric properties match those of conventional mask thermoplastics, ensuring consistent coverage and accurate dose delivery.



R-601 1lb jar	R-603 3lb jar
-------------------------	-------------------------

ThermoSheets™

Klarity ThermoSheets™ are thermoplastic sheets for bolus applications, available in both solid and perforated patterns. They can be molded like a standard thermoplastic mask, making them ideal for treating irregular contours. Heated in an oven or water bath, ThermoSheets™ are easy to handle and form, with dosimetric properties matching conventional thermoplastic masks.



Solid	
R-1600A 46cm x 61cm 1.6mm	R-4000A 46cm x 61cm 4.0mm
R-3200A 46cm x 61cm 3.2mm	R-4000E 15cm x15cm 4.0mm
R-3200 92cm x 61cm 3.2mm	R-4000F 30cm x 30cm 4.0mm
Perforated	
R-2436A 46cm x 61cm 36% perforation 2.4mm	R-3242A 46cm x 61cm 42% perforation 3.2mm
R-2436 92cm x 61cm 36% perforation 2.4mm	R-3242 92cm x 61cm 42% perforation 3.2mm

PATIENT SPECIALIZED ACCESSORIES

SHIELDS

Eye Shields

Klarity offers Eye Shields with lower transmission than other eye shields to protect the patient's lens during radiation therapy. Each Tungsten Eye Shield comes with both a 0.5 mm and 1 mm anodized aluminum cap to reduce electron backscatter. The shield can also be used without a cap for more superficial treatments. Please visit our website for full product dimensions and specifications.

KRP-936-583 Extra-Small Tungsten Eye Shield with Aluminum Caps 2mm	KRP-936-598 Small Tungsten Eye Shield with Aluminum Caps 3mm	KRP-936-589 Large Tungsten Eye Shield with Aluminum Caps 2mm	 KRP-936-627 Extra-Large Tungsten Eye Shield with Aluminum Caps 3mm
KRP-936-596 Extra-Small Tungsten Eye Shield with Aluminum Caps 3mm	KRP-936-587 Medium Tungsten Eye Shield with Aluminum Caps 2mm	KRP-936-623 Large Tungsten Eye Shield with Aluminum Caps 3mm	
KRP-936-585 Small Tungsten Eye Shield with Aluminum Caps 2mm	KRP-936-601 Medium Tungsten Eye Shield with Aluminum Caps 3mm	KRP-936-591 Extra-Large Tungsten Eye Shield with Aluminum Caps 2mm	

Testicle Shield Stand

Both sets of shields should be used with the adjustable stand.



KRP-926-100
Adjustable Testicle Shield Stand
25.4cm x 12.7cm x 34.29cm
10.6lbs

Lead Testicle Shields

The lead testicle shields are designed with a tongue and groove interlock to eliminate radiation leakage. The top and bottom halves can be secured together with two rubber straps (supplied with each unit) or with nylon tape. Rotating the entrance of the shield toward the ceiling will reduce scatter radiation from entering the shield. No rubber straps are needed if using the adjustable stand.



KRP-928-200
Small Lead Testicle Shield
Inner Diameter: 5.09cm
Outer Diameter: 7.65cm
Opening: 1.53cm x 3.19cm
1.27cm lead wall

KRP-928-250
Medium Lead Testicle Shield
Inner Diameter: 6.37cm
Outer Diameter: 8.92cm
Opening: 1.91cm x 4.08cm
1.27cm lead wall

KRP-928-300
Large Lead Testicle Shield
Inner Diameter: 7.65cm
Outer Diameter: 10.2cm
Opening: 2.55cm x 5.09cm
1.27cm lead wall

CT Testicle Shields

CT Testicle Shields are made of blue silicone rubber with shore hardness of 50A.



KRP-926-200
Small CT Testicle Shield
Inner Diameter: 5.09cm
Outer Diameter: 7.65cm
Opening: 1.53cm x 3.19 cm

KRP-926-250
Medium CT Testicle Shield
Inner Diameter: 6.37cm
Outer Diameter: 8.92cm
Opening: 1.91cm x 4.08cm

KRP-926-300
Large CT Testicle Shield
Inner Diameter: 7.65cm
Outer Diameter: 10.2cm
Opening: 2.55cm x 5.09cm

PATIENT SPECIALIZED ACCESSORIES

SKIN MARKING

Steritatt

Klarity offers SteriTatt, the first medical marking device created specifically for radiation therapy. This single-use solution supports precision treatment and features sterile, non-toxic, low-allergy ink.


KST-90001 SteriTatt Black Ink Cartridge 25 per box	KST-90015 SteriTatt Pink Ink Cartridge 25 per box
KST-90002 SteriTatt Black Ink Cartridge with Needle 25 per box	KST-90050 SteriTatt Black Ink Dropper. 3CC. 10 droppers per box



SureMark®

Precision Skin Marking Systems

Klarity offers a full line of SureMark® labels and markers, including nonmetallic, semi-lucent, and liquid fiducial markers. Please visit www.klaritymedical.com/suremark to see the full availability of SureMark® labels.

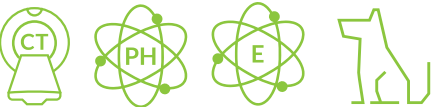




GENERAL TREATMENT ACCESSORIES

COMPRESSION BELT



R624-SCF1-5A
Compression Belt
183cm x 22cm

- Easy to clean
- Adjustable to accommodate many different body types
- Easy-to-use release button on the pump
- All-in-one design with no assembly required

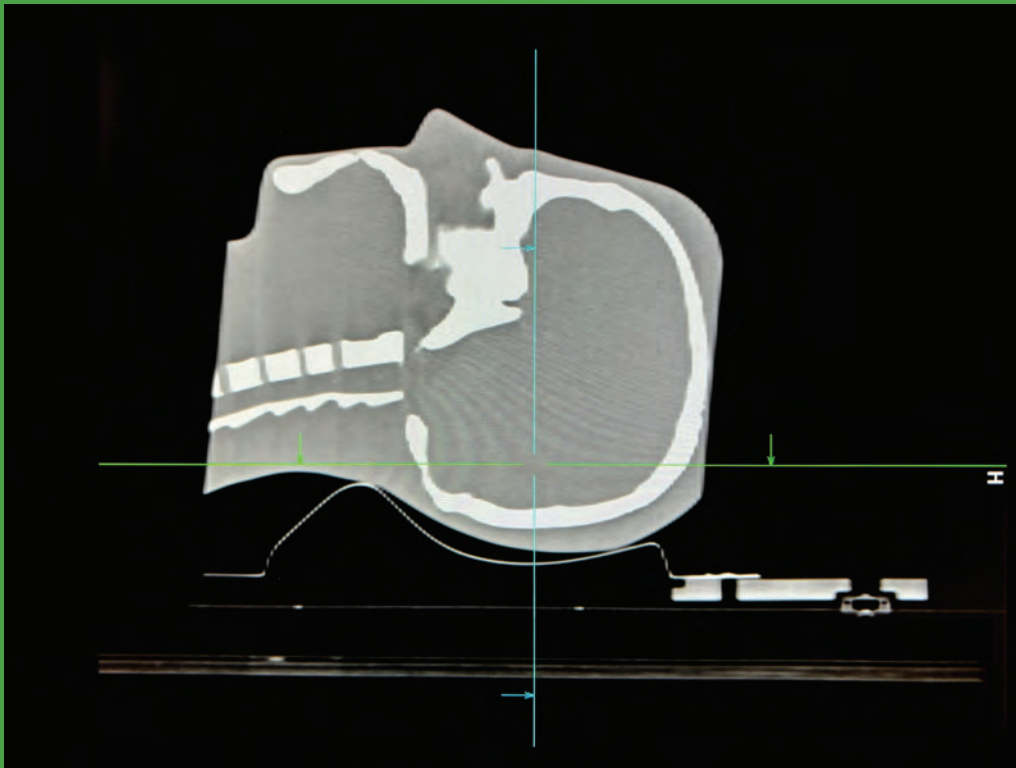

Carbon Fiber Head Support

GENERAL TREATMENT ACCESSORIES

HEAD SUPPORTS

Carbon Fiber Head Support

Klarity Carbon Fiber Head Supports deliver superior strength and a lightweight design, easily indexing with any standard Silverman pins. Order them individually or as a set of six to cover every height your setup requires.



- R507-CF**
Carbon Fiber Head Support Set of 6
- R507-ACF**
Carbon Fiber Head Support A
- R507-BCF**
Carbon Fiber Head Support B
- R507-CCF**
Carbon Fiber Head Support C
- R507-DCF**
Carbon Fiber Head Support D
- R507-ECF**
Carbon Fiber Head Support E
- R507-FCF**
Carbon Fiber Head Support F

Silverman Head Support

Klarity Clear Silverman Head Supports fit a variety of couches and tabletops. Order them individually or as a set of six for a complete range of height options.



- R507**
Silverman Head Support Set of 6
- R507-A**
Silverman Head Support A
- R507-B**
Silverman Head Support B
- R507-C**
Silverman Head Support C
- R507-D**
Silverman Head Support D
- R507-E**
Silverman Head Support E
- R507-F**
Silverman Head Support F

Prone Head Support

Klarity Prone Head Support provides comfortable prone positioning by supporting the forehead, cheeks, and chin while reducing pressure around the eyes, nose, and mouth. Its cushion is mounted on an acrylic plate that indexes to standard Silverman pins, offering compatibility with any baseplate using Silvermans.



R520
Prone Head Support



GENERAL TREATMENT ACCESSORIES

INDEXING BARS




RNE-100AL

Turn-Lock Aluminum Indexing Bar
3.2mm thickness
Compatibility: Exact™ Couch Top



RNE-102AL

Press-Lock Aluminum Indexing Bar
3.2mm thickness
Compatibility: Exact™ Couch Top

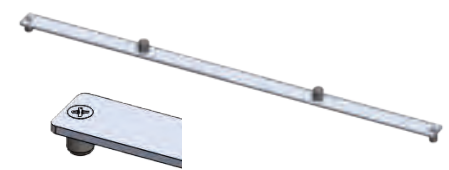

RNE-103AL

Insert-Lock Aluminum Indexing Bar
3.2mm thickness
Compatibility: iBEAM®, Novalis, and TomoTherapy® Couch Tops



RNE-104AL

C-Lock Aluminum Indexing Bar
3.2mm thickness
Compatibility: Civco GE CT Couch Top



RNE-104CLP

Low-Profile C-Lock Aluminum Indexing Bar
3.2mm thickness
Compatibility: Civco GE CT Low-Profile Thin Couch Top



RNE-106AL

Press-Lock Aluminum Indexing Bar
3.2mm thickness
Compatibility: Exact™ 48cm Couch Top



RE-6E

Universal Aluminum Indexing Bar
3.2mm thickness
490-535 mm range
Compatibility: All Couch Tops




RE-7

Aluminum Indexing Bar
6mm thickness
Compatibility: Varian Truebeam Couch Tops

Optek™ Indexing Bars

Compatible with the
Klarify Optek™ Overlay

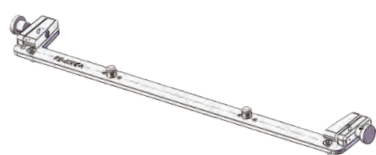


RE-630EA

Insert-Lock Aluminum Optek™ Indexing Bar
6mm thickness
Compatibility: iBEAM®, Novalis, and TomoTherapy® Couch Tops



RE-7EW

Turn-Lock Aluminum Optek™ Indexing Bar
6mm thickness
Compatibility: Varian Couch Tops

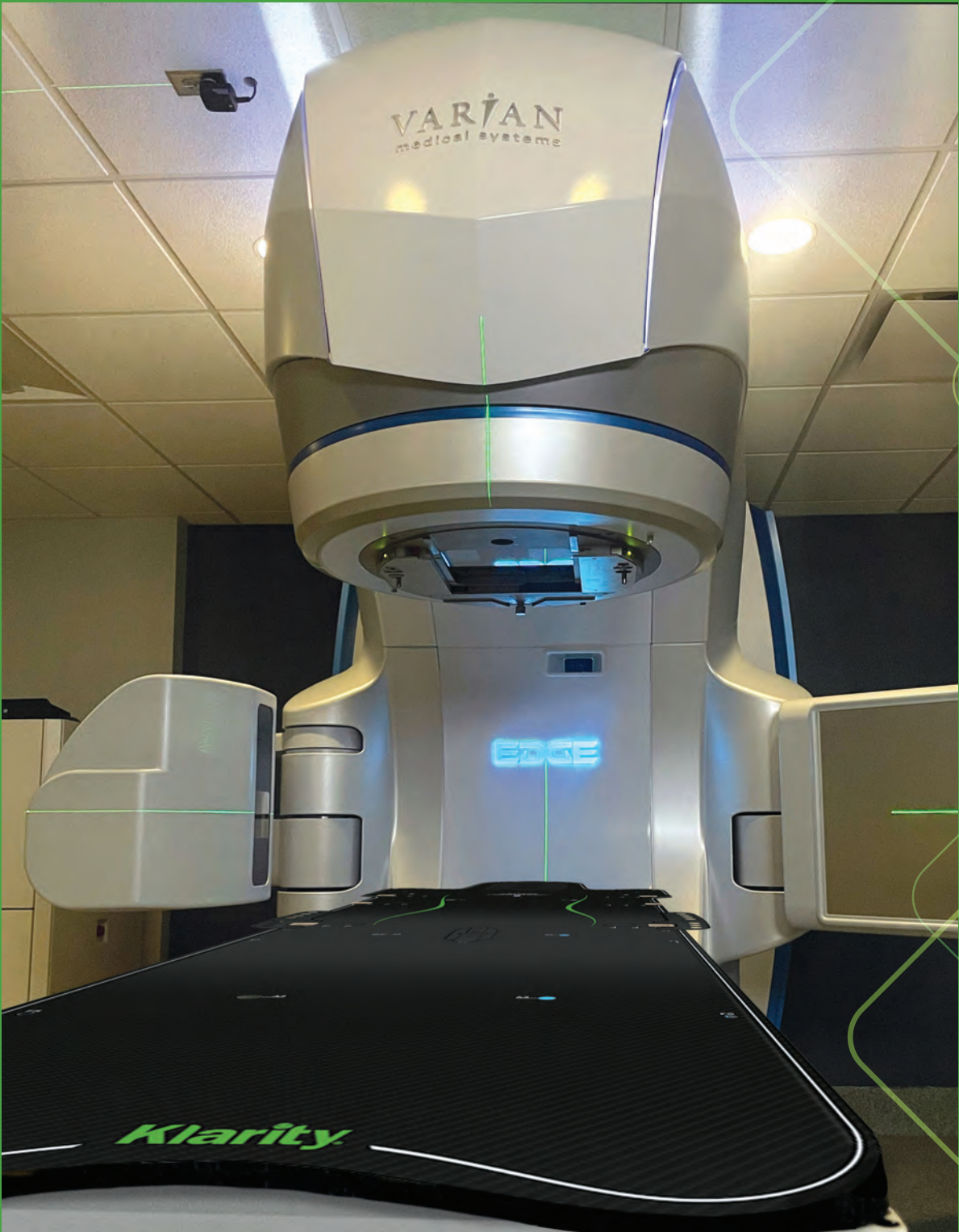

GENERAL TREATMENT ACCESSORIES

TABLE PADS & SUPPORT BLOCKS

Klarity's table pads, wedges, blocks and straps are engineered for long-term, multi-patient use. They feature durable, easily-cleanable surfaces that slide smoothly beneath patients while maintaining stability.

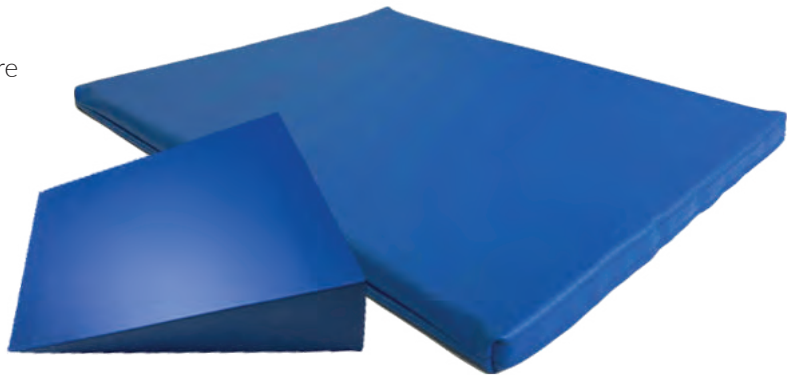


Table Pads



KRA1-5266
Short Table Top Pad
52cm x 66cm
25mm



KRA1-52100
Medium Table Top Pad
52cm x 100cm
25mm



KRA1-52150
Long Table Top Pad
52cm x 150cm
25mm

Support Blocks



KFB-005
Support Block 1
18cm x 25.5cm
2.5cm



KFB-006
Support Block 2
18cm x 25.5cm
5cm



GENERAL TREATMENT ACCESSORIES

WEDGE CUSHIONS & STRAPS


Wedge Cushions



KFB-001
Wedge Cushion 1
51cm (length) x 25.5cm (width)
x 18cm (height)
35° & 55°



KFB-002
Wedge Cushion 2
51cm (length) x 40.6cm (width)
x 30.5cm (height)
37° & 53°

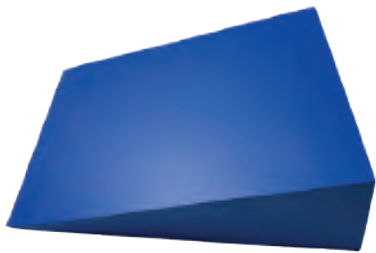

KFB-003
Wedge Cushion 3
51cm (length) x 76cm (width)
x 40.6cm (height)
28° & 62°



KFB-004
Wedge Cushion 4
25.5cm (length) x 18cm (width)
x 7.5cm (height)
17° & 73°



R516-229L
Green Bolster
17in (length) x 8in (width)
x 6.5in (height)

Straps



R565-603C
Green Cleanable Straps
183cm x 12.7cm



R565-607C
Black Cleanable Straps
183cm x 25.5cm

Supporting the Future of Cancer Care

At Klarity, we believe in investing in the next generation of Radiation Therapists. By donating equipment to schools and training programs, we give students the opportunity to practice, learn, and build confidence before stepping into real-world clinical settings.

“I want to express our sincere gratitude to you and the Klarity team for your incredible generosity. Your support means so much to our program, and we are genuinely excited to integrate these resources into our classroom and lab settings. These tools will undoubtedly enhance the educational experience for our students and help them prepare for success in the field.

We truly appreciate that Klarity has a program dedicated to supporting educational initiatives like ours. It speaks volumes about your commitment to the future of radiation therapy and the professionals who will drive it forward.”

—Angela Oliveira, MPA, RT(T)
Kakos School of Arts and Sciences
Riverdale, NY

Klarity is proud to support students as they develop the skills that will help improve patient care for years to come.

For more information, email us at
info@klaritymedical.com

OVENS

IN THIS SECTION:

Ovens..... 98

Oven Accessories.....100


KLARITY OVENS

A better, safer heating method

Dry-Heat Technology for Thermoplastics

Klarity pioneered dry-heat technology to efficiently prepare all standard thermoplastic masks, Klarity Cushions™, and bolus materials. By eliminating water, dry-heat ovens heat devices faster, reduce cross-contamination risk, and offer warm, dry devices that patients find more comfortable.

The Klarity AirFlow™ Oven

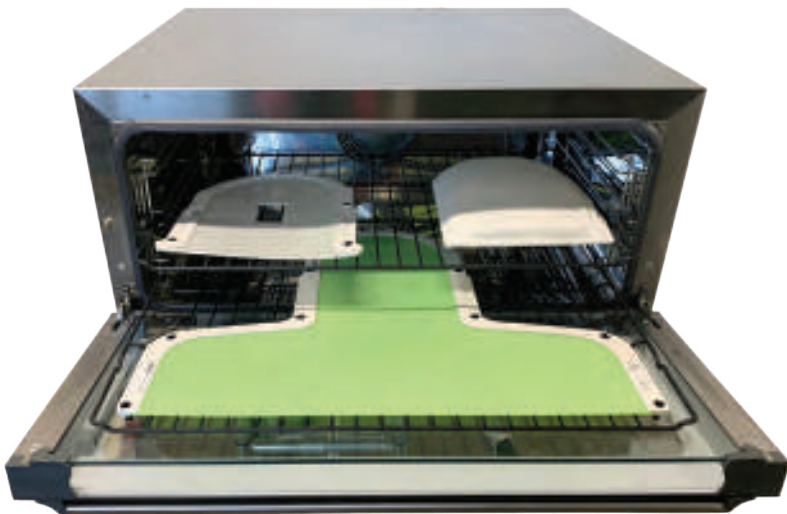

KEL-2400	KEL-1400
Klarity AirFlow™ Oven (220V)	Klarity AirFlow™ Oven (110V)
Outer dimensions: 88cm (length) x 80cm (width) x 41cm (height)	
Inner dimensions: 48cm (length) x 64cm (width) x 27cm (height)	
Weight: 40kg	



Why Choose Dry Heat?

- **Efficiency:** Ovens heat up quickly and eliminate the need for drying, streamlining workflow
- **Safety:** No hot water means reduced bacterial risk and lower chance of cross-contamination
- **Comfort:** A dry, warm mask is gentler on patients than a hot, wet mask

Faster Heating

Reaches 165°F in under ten minutes. Masks heat in just 3–5 minutes, while Klarity Cushions™ are ready in about 10–12 minutes.

Simple, Flexible Operation

A digital interface allows precise temperature control. The high-powered fan ensures consistent heat distribution.

Designed for Multiple Devices

Double racks accommodate Head and Shoulder Masks or multiple smaller devices. Non-stick Teflon coating means you can load masks directly with no fuss.


Easy to Clean and Move

A wipe-down interior and optional stainless-steel cart simplify maintenance and mobility.

“Once we received our Klarity oven, we put our water bath in storage and never used it again!”

—Sarah Steffen, BAS., R.T.(T)
Swedish Cancer Institute – Issaquah
Issaquah, WA

OVEN ACCESSORIES

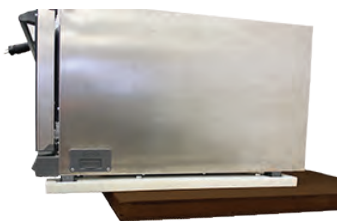

KEL-113

AirFlow™ CounterMount
Countertop extender for the AirFlow Oven
66cm (length) x 7.6cm (width) x 2.5cm (height)



KEL-1200RK

Supplemental Oven Rack



KEL-LINER

Non-Stick Teflon Oven Liners
43cm x 63.5cm



KEL-C100

AirFlow™ Oven Cart
91.5cm (width) x 76cm (depth) x 71cm (height)
Weight: 23kg

INCLUDED:

The CounterMount comes as a complete kit with two HDPE counter extenders, all hardware, multiple bolts for varying countertop thicknesses, and an installation template

Explore More Oven Resources

- **Scientific Article:** Learn about waterbath contamination in the 2014 *Radiation Therapist* article by Brewer, Ravine, and Bru.
- **Instructional Video:** Watch a step-by-step guide on preparing devices with dry-heat ovens.
- **CE Class:** Earn continuing education credits while mastering best practices in thermoplastic preparation.
- **FAQ:** Visit www.klaritymedical.com/ovens for answers to common questions.